Dreame Nova IPL Home Use Hair Removal Device

User Manual

DREAME

EN	User Manual.....	001
FR	Manuel d'instructions	023
ES	Manual de usuario	042

Warning

Please read this manual carefully.
Improper operation may cause
personal injury or property damage.

Warning sign description: Do not use this device while bathing or showering.

Do not use this device near bathtubs, sinks, or other water-filled containers. Store the device and power adapter in a dry place.

If the device is damaged, do not touch the internal components to prevent electric shock.

A damaged power cord must be replaced by the manufacturer, authorized service department, or qualified professionals to prevent danger.

Do not use unofficial adapters, and do not place this device improperly.

Do not use this device near the eyes to prevent eye injury.

Always wear goggles when operating the device.

For back hair removal, it is recommended to have someone assist you to avoid risks.

Adjust the energy intensity according to the skin tone and condition of each body part. Always begin with L1 for the first time.

No data is available for use on children under 15. Keep the device away from them. Operate at a safe distance from children.

Do not apply antiperspirant immediately after treating underarms.

Regularly check and clean the aperture to prevent skin burns or device malfunctions.

Keep the aperture away from objects that are easily demagnetized to prevent demagnetization or device malfunction.

Do not use the device if the vents are blocked to prevent fire.

Do not aim the flash at black objects to prevent device malfunction.

Do not use the device on clothing or areas not intended for treatment, to prevent clothing discoloration, material damage, or singeing.

Do not use the device on pets to prevent skin burns, singed fur, or stress responses.

After the device is powered on, if the built-in fan does not work (no fan operation sound), do not use the device to prevent smoking or fire.

Regularly clean residual hair and skincare product from the device to prevent device malfunctions or skin burns.

This appliance shall not be used by children and shall not be applied on the skin of children. Children shall not play with the appliance.

These appliances are not to be used on children; children being supervised not to play with the appliance.

This appliance can be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

For a part of class III construction supplied from a detachable power supply unit, the appliance is only to be used with the unit provided.

WARNING. Possible eye injury (potentially leading to loss of vision) or skin injury if instructions are not followed. Read and follow the instructions.

WARNING. Damage of the skin may occur after prolonged or repeated surface application on one site.

Misuse can lead to eye damage. Protect the eye from exposure. Do not use the appliance over the eyelids or close to the eye.

WARNING. Do not override the safety mechanisms inherent to the device.

WARNING. If taking a photo sensitizing medication or herbal remedy, seek medical advice prior to application of the appliance as use may harm the skin.

Do not use the appliance in areas where skin integrity has been impaired (for example, in the presence of burns, lesions, blisters, scars, cuts, open wounds, active skin disease, recent sun tan or sun burn, and/ or infection).

This appliance is not intended for the treatment of medical conditions. Medical advice should be sought to address, for example, moles, skin rash, itchy skin, skin fungus or infection, skin bumps, or skin tags.

For users with different skin tones and hair colors, please refer to the Skin Tone and Hair Color Table in this manual.

Avoid over-usage which includes excessive passes, stacking of pulses and/or increased frequency of use.

Verify that there is no visible damage to the emission optic or the appliance housing. If there is damage, do not use this appliance.

For instructions on how to test skin reaction prior to full-area treatment/use, please refer to the following chapter: Skin Test.

Before use, prepare the skin—for example, by cleaning hair and removing make-up, jewelry, and piercings.

For other relevant information, please refer to the following chapter: Instructions on Device Operation, Disinfection, Cleaning, Maintenance, and Storage.

Before use, please read the instruction manual carefully, then if you still have any questions or uncertainties, please consult a doctor for advice first.

Precautions


Individuals with naturally dark skin should not use this device

Avoid excessive sun exposure before hair removal

Shave the hair before using this device (do not use chemical hair removal methods)

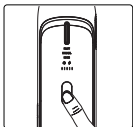
Always begin with L1 for the first use

Do not use this device in dark environments

Always wear goggles while using the device

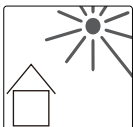Do not use this device under direct sunlight





Protect the device from severe impacts or shaking

Product Overview

This device transforms the standard experience with its built-in ice-cooling system, instantly soothing heat and minimizing skin irritation while removing hair. Its faster flash speed also makes hair removal quicker and easier.


Before using this device, shave the target area (do not use chemical hair removal methods). The device uses intense pulsed light (IPL) that goes into the skin. The light gently heats the follicle to change its structure, stopping hair growth at the source. This leads to hair naturally falling out for long-lasting results.

Components

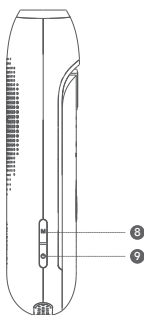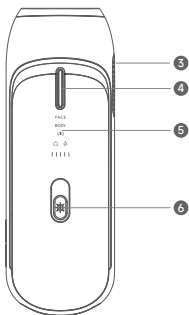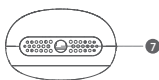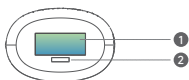
- 1 Aperture
- 2 Skin tone sensor indicator
- 3 Air outlet
- 4 Sensor indicator
- 5 Mode display
- 6 Flash button
- 7 DC power port
- 8 Mode switch button
- 9 Power button/intensity switch button
(power/intensity button)

Accessories are shown below: adapter, razor, goggles, storage bag

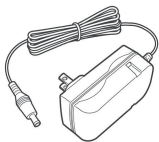Adapter

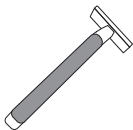Razor


Goggles

Storage bag

Indicators and Modes

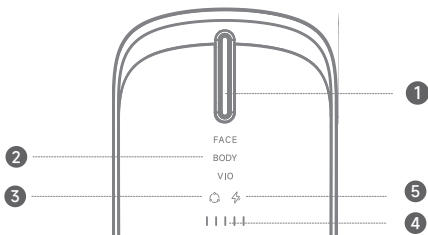- 1 " | " Sensor indicator: This remains off when the flash conditions are not met, and stays on when the conditions are met.
- 2 " FACE " Face mode indicator: This is the default mode displayed when the device is turned on.
" BODY " Body mode indicator.
" VIO " Bikini mode indicator.
- 3 " ⌚ " Smart skin tone detection indicator: Smart skin tone detection is enabled by default after the device is turned on. When the aperture is placed vertically against the skin, the device will automatically determine the appropriate light intensity through the skin tone detection feature, and the intensity indicator will illuminate. When the device is lifted away from the skin, all intensity indicators will turn off.
- 4 " ||||| " Intensity indicator: 5 intensity levels are available. After the device is turned on, press the "power/intensity button" to cycle through automatic skin tone detection mode and manual intensity settings in the following order: automatic skin tone detection - 1 - 2 - 3 - 4 - 5 - automatic skin tone detection...
- 5 " ⚡ " Continuous flash mode indicator. The device is set to manual flash mode by default. Press and hold the "flash button" for 2 seconds to switch to continuous flash mode, and this indicator will illuminate. In continuous flash mode, press and hold the "flash button" for 2 seconds again to switch to single flash mode.

Note: The device will automatically switch to manual flash mode if no skin contact is detected for 30 seconds.

Instructions for Use

Before using the device

Thoroughly clean the treatment area. Make sure the skin is free of hair and oil residue.

Check the device's aperture. Make sure it is free from dust, hair, or other debris, and has no damage or cracks.

Intensity level settings

This product is designed for human hair removal with 5 available intensity levels. You can refer to the skin tone comparison card to select the appropriate intensity level.

Note:

1. Place the skin tone comparison card against your skin and select the appropriate intensity level based on your skin tone and hair color.
2. If you feel pain or discomfort during use, lower the intensity level. Individuals with darker skin tones absorb more light and are advised to use lower levels. During normal use, the light should only produce a mild warming sensation on the skin, not pain.
3. Please keep this card in a safe place for future reference.

Skin tone and hair color table

Skin tone / Hair color		Skin tone					
		White	Light	Medium	Dark	Very Dark	Black
White/gray		X	X	X	X	X	X
Light blond		X	X	X	X	X	X
Red		X	X	X	X	X	X
Dark blond/ light brown		✓	✓	✓	✓	X	X
Brown		✓	✓	✓	✓	X	X
Dark brown		✓	✓	✓	✓	X	X
Black		✓	✓	✓	✓	X	X

(Note: ✓ means suitable, X means not suitable)

Note: Make sure your skin type is suitable for using this device, and refer to the skin tone comparison card for verification. Skin reactions may vary when treatments are performed at different times and in different settings.

Skin test

If this is your first time using the device, or if you have a recent tan, please perform a skin test first. This checks your skin's reaction and helps set the right intensity level.

1. Select a small area within the treatment zone.

Note: Do not test on sensitive or bony areas such as ankles.

2. Press and hold the power button for 1 second and select the L1 intensity level.

3. Place the device vertically against the skin, with the aperture pressed onto the skin.

4. Press the flash button to emit a flash.

5. Keep the device in contact with the skin, and move to the next test area.

6. Gradually increase the intensity level, allowing the device to flash once each time you change to a higher level, then move to the next area. Repeat this operation for all recommended intensity levels for your skin type.

7. Observe the skin condition two hours after the skin test. If no abnormal symptoms such as redness, swelling, or blistering are observed in the test area, you can proceed to use this device.

Hair removal guide for different areas

1. Underarms

- Before hair removal, clean the underarms to remove any residue of antiperspirants, perfumes, creams, or other cosmetic products.
- The skin in the underarm area can be loose. For areas with loose skin, stretch the skin with your fingers before hair removal to keep the skin flat.
- A slight heating sensation is normal during flashing.

2. Bikini line

- The skin in this area is sensitive. For darker or highly pigmented areas, start from L1 and allow the skin to adapt gradually.
- For areas with loose skin, stretch the skin with your fingers before hair removal to keep the skin flat.
- A slight heating sensation is normal during flashing.

3. Hands and feet

- For hands and feet, where bones are closer to the surface and the skin may be uneven (e.g., around joints), you might feel more pain during treatment.

Before Using the device

1. Pre-Use Skin Check

Ensure the treatment area has no broken skin, severe sunburn, redness, or rashes.

Avoid if you've recently done high-irritation procedures (like laser or chemical peels) on the same area.


Consult your doctor if pregnant, photosensitive, or using photosensitizing medications.

2. Shave First

Use a razor to shave the target area down to skin level.

Avoid waxing or plucking, as this can prevent light energy from reaching the hair root effectively.

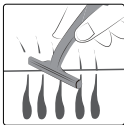3. Cleaning & Dry

Clean the flash area with clean water or mild cleaning products.


Dry completely to keep it dry, oil-free and unobstructed.

4. Connect to Power

Plug the adapter into a power outlet and connect it to the device's DC port.

Place in a ventilated, dry area away from water sources.

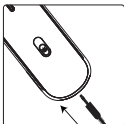Using the device

1. Power On

Press and hold the Power/Intensity Button for about 1 second to turn on the device.

Default settings after power-on:


- FACE Mode
- Auto Skin Tone Detection Mode ( lights up)
- Cooling function activated

2. Choose Your Mode (FACE / BODY / VIO)

Press the Mode Button briefly to switch between the three modes below:

- FACE Mode: Suitable for soft and fine hair areas such as the face, lip line and fingers
- BODY Mode: Suitable for large body areas such as the arms, calves and underarms
- VIO Mode: Suitable for private and sensitive areas such as the bikini line

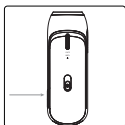Note: Start with FACE/BODY Mode first. Try the VIO area after getting used to the device.

3. Pre-cool

The light-emitting area cools down automatically after power-on. Wait 10–20 seconds, then gently touch it with the back of your hand to confirm it feels noticeably cool. This ensures a more comfortable treatment.

US 4. Wear Goggles

NEVER look directly at the flash during use. Goggles must be worn at all times.

5. Skin Tone Recognition ()

The device powers on in this mode automatically. () will light up, no intensity level is displayed.

Press the treatment window firmly against skin. After successful skin tone recognition, the recommended intensity level (1-5) appears automatically. Press the flash button to treat.

Lift the device off your skin, only () lights up, indicating it remains in Skin Tone Recognition Mode.

6. Choose Your Intensity

Press the Power/Intensity Button once to switch from Auto Mode to Manual Mode. Your current level (1-5) will be displayed.

Levels will cycle in this order: Auto → 1 → 2 → 3 → 4 → 5 → Auto →

...

It is recommended to start at a low level, then increase it step by step based on your own comfort level.


To go back to Automatic Skin Tone Detection Mode, keep pressing the button briefly until the device switches back to Auto ().

7. Manual Flash & Continuous Flash


Manual Flash Mode

- Hold the treatment window vertically and firmly against the target hair removal area — the sensor light strip will light up.
- Press the flash button briefly once, and the device will complete a full round of flashes in the current position.
- FACE / VIO: Each press triggers 2 consecutive flashes. Wait until the flashing stops before moving to the next area.
- BODY: Each press triggers 4 consecutive flashes.

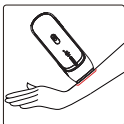Continuous Flash Mode

- Hold the treatment window vertically and firmly against the target hair removal area. Press and hold the flash button for about 2 seconds — the ⚡ will light up, indicating that Continuous Flash Mode is activated.
 - FACE / VIO: Each press → device automatically delivers 2 consecutive flashes.
 - BODY: Each press → device automatically delivers 4 consecutive flashes.
- Press and hold the flash button for about 2 seconds again to exit Continuous Flash Mode.
- If no skin contact is detected within 30 seconds, the device will automatically exit this mode.

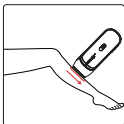Contraindications

1. Pregnant women, breastfeeding women, women during menstruation, and individuals with scars on the skin, diabetes, sensitive skin, or photosensitivity should not use this device.
2. Individuals with severe heart disease, severe hypertension, malignant tumors, folliculitis, infectious skin diseases, skin burns, or blood coagulation disorders should not use this device.
3. Individuals with a history of collagen metabolism disorders (including scarring and poor wound healing) should not use this device.
4. Individuals should not use this device on areas with vascular diseases such as varicose veins and telangiectasia.
5. Individuals with photosensitive skin or are prone to skin rashes or skin discomfort should not use this device.
6. Individuals should not use this device on areas with infection, herpes simplex, skin damage, or hematoma.

7. Individuals with immune suppression-related diseases (including HIV infection, lupus erythematosus, and porphyria) should not use this device.
8. Individuals with bleeding disorders, or those using topical isotretinoin, alpha hydroxy acids (AHA), beta hydroxy acids (BHA), or azelaic acid should not use this device.
9. Individuals who have used chemical hair removal methods within 30 days, or have recently used whitening products or cosmetics containing vitamin A should use this device with caution.
10. Do not use this device with any volatile essential oils. Make sure that the skin in the treatment area is clean, dry, and free from creams, lubricants, or other cosmetics. If the skin is not clean and dry, it may affect the product efficacy or even damage the device.
11. Individuals taking medications such as aspirin, immunosuppressants, or isotretinoin should not use this device. Individuals who have recently taken photosensitive medications such as tetracycline or sulfonamides should not use this device.
12. Individuals with sunburned skin should wait until the damaged skin has fully recovered before using this device. This product should not be used on sensitive areas such as mucous membranes.
13. Do not use this device on skin with any unhealed wounds. This device is not recommended for individuals whose excessive hair growth is caused by hirsutism or other medical conditions.
14. Painkillers can make your skin less sensitive to heat. Use this device with caution before the medication wears off.
15. If your skin or hair color is not suitable for this product, please do not use it. (Consult a doctor before use if needed.)

Prohibited areas

1. Inside the mouth, private parts, mucous membranes	8. Eyebrows and hair
2. Areas that have undergone plastic surgery or other surgical procedures	9. Ears
3. Eyeballs, eyelids, and areas around the eyes	10. Nipples and areolae
4. Adam's apple	11. Inside the navel
5. Head	12. Areas with large veins
6. Areas with acne or pimples	13. Areas with tattoos or permanent makeup
7. Areas with pigmentation or moles	

Note: This list is not exhaustive. Consult a doctor if you are uncertain about treating any specific area.

Care

Treatment frequency

Weeks 1-4: 2-3 treatments per week.

Weeks 5-8: 1-2 treatments per week.

After week 8: 1 treatment per month or 1 treatment every two weeks.

Note: Adjust the treatment frequency based on your hair growth. For thicker hair, you may need to increase the intensity level and treatment frequency. This varies by individuals.

After using the device

1. Avoid contact with water within 6 hours after hair removal. After 6 hours, you may use toner to add moisture to the treatment area.
2. After hair removal, do not use perfume, face masks, or cosmetics containing exfoliating particles on the treatment area.
3. Avoid UV exposure within 48 hours after hair removal.
4. Sunless tanning is prohibited during the course of your hair removal treatments.
5. If you plan to undergo plastic surgery on the treatment areas after hair removal, please consult your doctor.

Adverse reactions and handling methods


Adverse reaction	Description	Handling method
Dry or itchy skin	The hair removal process slightly removes oil and moisture from the skin surface.	Discomfort will disappear in 1 day. For quick relief, you can apply a mild, non-irritating moisturizer.
Pain	The level of pain varies from person to person and is related to skin tolerance, the intensity level used, and the treatment area (e.g., bony spots are more sensitive).	Lower the intensity by 1-2 levels and try again. If pain persists, it is recommended to stop using the device and wait for the skin to recover before trying again.
Mild redness or burning sensation	This is a normal skin reaction after light exposure.	No special treatment is needed. The redness or burning sensation will disappear over time. Alternatively, you can apply a cooling gel or a cold compress.
Burnt hair	Hair is not thoroughly shaved before hair removal. Residual hair heats up under the light, causing a slight burning smell.	Before the next use, thoroughly shave the treatment area with a razor.

Cleaning and Maintenance

To ensure optimal product performance and extend its service life, clean the device before and after each hair removal session to prevent accumulation of hair and dust. Lack of regular cleaning may affect hair removal effectiveness and even shorten the device lifespan.

Check the aperture and device components to ensure they are clean. Do not rinse the device and accessories under running water. Do not use steel wool, abrasive cleaning agents (such as toothpaste), or corrosive liquids (such as gasoline, engine oil, and acetone) to clean the device.

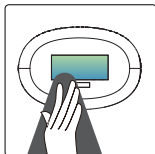Press and hold the Power/Intensity Button for about 2 seconds to turn off the device. Gently wipe the flash window and device body with a dry, lint-free soft cloth. Place the device and adapter back into the original packaging. Store in a well-ventilated, dry area, away from direct sunlight and out of reach of children.

Note: Auto Power-Off activates after 5 minutes without operation.

Product Specifications

Product name	Hair Removal Device
Model	R3C16-G Pro
Product dimensions	7.28x2.56x1.77inches
Net weight	About 0.705pounds
Package contents	Hair removal device × 1, user manual × 1, adapter × 1, Goggles × 1, razor × 1, storage bag × 1
Adapter's rated input	100-240V~ 50Hz/60Hz,1.5A Max
Output voltage/current	24 V  3 A
Rated power	72 W
Operation environment temperature	41°F~86°F (5°C~30°C)
Intensity levels	5 levels
Main filter wavelength	510 nm-1200 nm
Place of manufacture	China

Possible Side-Effects

There is few side-effect when using Hair Removal Device according to the instructions for use. However, every beauty process, including the usage of home-use devices, may cause some degree of side-effects. If any of the following conditions (refer to Side Effect Table below) appear and last more than 24 hours, it is suggested to consult your physician or dermatologist.

Side Effect Table:

Degree of Side Effects	Side Effects
Minor Side-Effects	Irritation Redness Sun Sensitivity
Rare Side-Effects	Burns Scarring Pigment Changes Excessive Redness Swelling
Extremely Rare Side-Effect	Infection Bruising

Maintenance and Storage

Cleaning

- Make sure to turn off and unplug Hair Removal Device before cleaning.
 - It is recommended to clean the treatment window before and after every session by using a dry and clean cloth to gently wipe.
NOTE: If the dirt is difficult to be cleaned, you can dip a cotton swab or clean cloth in 75% medicinal alcohol to clean.
 - Keep Hair Removal Device dry. Never submerge the device or any parts of it in water, as it may lead to an electrocution risk.
 - Do not use steel wool, abrasive cleaners (toothpaste, etc.) or corrosive liquids (gasoline, motor oil, acetone, etc.) to clean the product.
- NOTE: When cleaning is finished, please conduct a visual inspection to check whether there is residual dirt. If the device is not visually clean, please repeat the above cleaning steps.

Storage

- To avoid breakage or missing device and its accessories, please store Hair Removal Device and its accessories in the gift box after each use.
- Clean the device before you store it and store in a safe and protective place, otherwise, it may affect the product function.

Trouble Shooting

If you experience any problem with Hair Removal Device, please refer to information below to try to resolve it:

Problem	Solution
The device cannot be turned on	Make sure that the power adapter is plugged properly into the device and the electrical socket.
The device cannot emit a light	<ol style="list-style-type: none"> 1) Verify whether the treatment window is in full contact with the skin. The flash button will be activated only when the device is in full contact with the skin 2) When the lamp tube reached its end, the device will not shoot a light.
If you feel overheat and pain during hair removal	<ol style="list-style-type: none"> 1) Shave and clean the hair in the expected treatment area before using. 2) Choose one comfortable energy level for yourself. 3) DO NOT treat the same area repeatedly at the same time.

If you still unable to solve a problem, please contact after-sales service.

Disposal

DO NOT throw away the Hair Removal Device with normal household waste at the end of its life, but hand it in to an official collection point for recycling (Contact your local town or city officials for recycling information). By doing this, you help to preserve the environment.

Warranty

The warranty period of the Hair Removal Device is 2 years from the date of purchase. Please contact local agency or after-sale service personnel if there is any problem during the period of warranty, repair according to relevant warranty clause.

If any operator requests more information such as circuit diagrams, parts list and product descriptions, for repairs carried out by qualified technical personnel, please contact us.

Symbols

Symbol	Meaning
	Batch code
	Manufacturer
	Date of manufacture
	Type BF applied part
	General warning sign
	"WEEE (Waste Electrical and Electronic Equipment)". The waste products should be handled legally.
	Refer to instruction manual/ booklet
	This way up
	Keep dry
	Fragile, handle with care
	Optical radiation
	Please wear protective glasses

Annex I: Manufacturer's EMC Statement

The device has no essential performance.

TABLE 1

Guidance and Manufacturer's Declaration – Electromagnetic Emissions	
Emissions test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class B
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies

TABLE 2

Guidance and manufacturer's declaration - electromagnetic Immunity		
Immunity Test	IEC 60601-1-2 Test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8kV contact ±2kV, ±4kV, ±8kV, ±15kV air	±8kV contact ±2kV, ±4kV, ±8kV, ±15kV air
Electrical fast transient/burst IEC 61000-4-4	±2kV for a.c. power lines	±2kV for a.c. power lines
Surge IEC 61000-4-5	±0.5kV and ±1kV for a.c. power line(s) to line(s)	±0.5kV and ±1kV for a.c. power line(s) to line(s)
Voltage dips, short interruptions and variations IEC 61000-4-11	0 % UT; 0.5 cycle at 0°,45°,90°, 135°, 180°, 225°, 270°, 315° 0 % UT ; 1 cycle 70 % UT; 25/30 cycle 0% UT; 250/300 cycle	0 % UT; 0.5 cycle at 0°,45°,90°, 135°, 180°, 225°, 270°, 315° 0 % UT ; 1 cycle 70 % UT; 25/30 cycle 0% UT; 250/300 cycle
Power frequency magnetic field IEC 61000-4-8	30 A/m 50Hz/60Hz	30 A/m 50Hz/60Hz
Conducted RF IEC61000-4-6	3Vrms 0,15 MHz – 80 MHz 6Vrms in the ISM frequency band	3Vrms 0,15 MHz – 80 MHz 6Vrms in the ISM frequency band
Radiated RF IEC61000-4-3	10 V/m 80 MHz – 2,7 GHz	10 V/m 80 MHz – 2,7 GHz
NOTE UT is the a.c. mains voltage prior to application of the test level.		

TABLE 3

Guidance and Manufacturer's Declaration - IMMUNITY to proximity fields from RF wireless communications equipment					
Immunity test	IEC60601 test level				Compliance level
	Test frequency	Modulation	Maximum power	Immunity level	
Radiated RF IEC 61000-4-3	385 MHz	**Pulse Modulation: 18Hz	1.8W	27 V/m	27 V/m
	450 MHz	*FM+ 5Hz deviation: 1kHz sine	2 W	28 V/m	28 V/m
	710 MHz 745 MHz 780 MHz	**Pulse Modulation: 217Hz	0.2 W	9 V/m	9 V/m
	810 MHz 870 MHz 930 MHz	**Pulse Modulation: 18Hz	2 W	28 V/m	28 V/m
	1720 MHz 1845 MHz 1970 MHz	**Pulse Modulation: 217Hz	2 W	28 V/m	28 V/m
	2450 MHz	**Pulse Modulation: 217Hz	2 W	28 V/m	28 V/m
	5240 MHz 5500 MHz 5785 MHz	**Pulse Modulation: 217Hz	0.2 W	9 V/m	9 V/m

Note* - As an alternative to FM modulation, 50 % pulse modulation at 18 Hz may be used because while it does not represent actual modulation, it would be worst case.

Note** - The carrier shall be modulated using a 50 % duty cycle square wave signal.

We, Dreame Trading (Tianjin) Co., Ltd., here by, declare that this equipment is in compliance with the applicable Directives and European Norms, and amendments. The full text of the EU declaration of conformity is available at the following internet address:
<https://global.dreametech.com>

WEEE Disposal and Recycling Information

All products bearing this symbol are waste electrical and electronic equipment (WEEE as in directive 2012/19/EU) which should not be mixed with unsorted household waste. Instead, you should protect human health and the environment by handing over your waste equipment to a designated collection point for the recycling of waste electrical and electronic equipment, appointed by the government or local authorities. Correct disposal and recycling will help prevent potential negative consequences to the environment and human health. Please contact the installer or local authorities for more information about the location as well as terms and conditions of such collection points.

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Manufacturer: Dreame Trading (Tianjin) Co., Ltd.