Dreame Astra IPL
Hair Removal Device
User Manual

| | | |
|-----------|----------------------------|-----|
| EN | User Manual..... | 001 |
| ES | Manual de usuario | 025 |
| FR | Manuel d'utilisation | 049 |

User Manual

IPL Hair Removal Device

Applicable Model: T32A

The IPL Hair Removal Device is an over-the-counter device indicated for the removal of unwanted hair. The device is also indicated for the permanent reduction in hair regrowth, defined as the long term, stable reduction in the amount of hair regrowing when measured at 6, 9 and 12 months after the completion of a treatment regime.

0 Welcome

Thanks for choosing IPL Hair Removal Device manufactured by Dreame Trading (Tianjin) Co., Ltd. In order to get better effect for IPL Hair Removal Device, we suggest you:

- Please read through the User Manual before using and keep this User Manual.
- For any questions or issues, please contact us at:
<https://global.dreame.tech.com/pages/contact-us>

Notice:

- This user manual guide you how to use IPL Hair Removal Device.
- We follow sustainable development strategy. So we keep the right that we will not notify customers when we change product function, interface and appearance.
- If there is difference between description and actual product due to upgrade, please refer to actual product.
- Please use device according to correct operation steps and treatment setting.
- Illustrations in this manual are for reference only.

Software release version and date: V1.1, February 09, 2026
User manual release version and date: V1.1, February 10, 2026

1. Safety Precautions

EN

Please first read these safety precautions carefully and then follow the user manual to operate. These safety precautions are to ensure the operation of product is correct and safe, to prevent users or other people from personal injuries or property damage, and please make sure to follow them.

1.1 Warnings

1)The IPL hair removal device shall not be serviced or maintained while in use with a PATIENT. The PATIENT is an intended LAY OPERATOR.

2)The isolated power supply is considered as part of the device.

3)Hair removal by lasers or intense pulse light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern, and South Asian heritage treated on the face and neck. This group should be careful.

4)DO NOT expose treated areas to the sun. Wait at least 24 hours after treatment before you expose treated skin to the sun or tan with artificial light. Even after 24 hours, make sure that the treated skin is not red from the treatment before you expose it to the sun. Be sure to carefully protect the treated skin with sunscreen throughout the hair removal period. Failure to follow this warning could result in skin swelling, blisters, or permanent skin discoloration or skin whitening.

5)DO NOT continue with the treatment if the treated area feels very hot, very sensitive or burns. Such use can result in burns or other skin injury. Stop treatment and call Contact Customer Service.

6)DO NOT use if you have received radiation therapy or chemotherapy treatments within the past 3 months.

7)DO NOT use if you have a history of skin cancer or areas of potential skin malignancies.

8)DO NOT use if you have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.

9)DO NOT use if you have an active implant, such as a pacemaker, incontinence device, insulin pump, etc.

10)DO NOT use if you have a history of herpes outbreaks in the area of treatment.

11)If your skin blisters or burns, STOP USE IMMEDIATELY! Continued use can result in burns or other skin injury.

12)The device should not be used on the face above the cheek line: around the eyes, eyebrows, or eyelashes, as this could cause serious eye damage.

13)Since the device is powered by the external power adapter, do not use the device by placing it in a location away from the electrical outlet.

14)Do not use the device together with electronic devices such as cell phone, electrotome, MRI device. Otherwise, it may affect the performance of the device or cause electric shock to the patient.

15)The applied part (AP) of the device is treatment window. In the process of use, the temperature of the AP may reach up to 60 °C. The operator can control the temperature by limiting the contact time, which is recommended to be less than one minute per session, and the operator can control the temperature limit by the contact time.

16)Keep out of the reach of children under 18 years of age.

17)If the power cable of the adapter becomes tangled, it may lead to strangulation, especially with excessive cable length.

18)Do not treat exactly the same area multiple times during one session to avoid skin redness.

19)Do not leave the equipment unattended when it is switched on to avoid the risk of fire or burns.

20)This equipment is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the equipment to avoid the risk of fire or burns.

21)This equipment is not intended for use by children. Children should be supervised to ensure that they do not play with the equipment to avoid the risk of fire and burns.

22)Men are prohibited from using this device for facial hair removal, as this may result in safety risks such as burns.

23)If the mains cord is damaged, you shall have it replaced by a service Centre authorized by Dreame Trading (Tianjin) Co., Ltd. or similarly qualified persons in order to avoid risk of electrocution.

24)Manufacturer will provide circuit diagrams, component part lists, descriptions, calibration instructions to assist to service personnel in parts repair.

25)Water and electricity are a dangerous combination! To avoid risk of electrocution,

- Do not use this equipment in wet surroundings (e.g. in the bathroom or near a shower or swimming pool).

- Do not let water run into the appliance.

26)Check if the voltage indicated on the equipment corresponds to the local mains voltage before you connect the equipment to avoid risk of electrocution or permanent damage to the equipment.

27)If the wall socket used to power the equipment has poor connections, the plug of the equipment becomes hot. Make sure you plug the equipment into a properly installed wall socket to avoid the risk of fire and burns.

28)The lay operator or lay responsible organization should contact the manufacturer or manufacturer's representative:

- for assistance, if needed, in setting up, using or maintaining the IPL hair removal device, or

- to report unexpected operation or events.

29)Keep easy to isolate from the supply mains

-Do not to position the equipment so that it is difficult to operate the disconnection device (plug).

30)WARNING:

-No modification of this equipment is allowed.

-Do not modify this equipment without authorization of the manufacturer.

-If this equipment is modified, appropriate inspection and testing must be conducted to ensure continued safe use of the equipment.

31)Pay attention to strangulation due to cables with excessive length, please be sure the cables are far away from neck when using the equipment.

1.2 Precautions

1)It is recommended to use this device in a well-lit indoor environment to avoid affecting the operation of the device.

2)When the IPL hair removal device is removed from the maximum or the minimum storage temperature environment to a room temperature environment (20°C), the IPL hair removal device is expected to take 2 hours to be ready.

3)When the treatment window of the device has been closely and correctly fitted to the expected skin area and the device fails to respond accordingly, the operator needs to restart the device. If the same problem still occurs after restarting, the built-in skin sensor may be loose or damaged; in this case, please contact our customer service.

4)During the operation, do not shield the air intake and ventilation holes to impact heat dissipation;

5)Do not throw away the device with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment;

6)Do not subject the equipment to heavy shocks to avoid risk of damage to the lamp;

7)Always unplug the appliance after use and in case of a power failure to avoid risk of damage to the equipment;

8)Keep the skin of the treated area clean, dry and free of any cosmetic residues before each operation.

9)During the treatment, try closing the device to restart if the display is not on. If the display is still not on after restart, please contact our customer service.

10)During the treatment, the trigger button can only be triggered when the treatment window is in good and correct contact with your skin. If the above specified conditions are not met and the trigger button can be triggered, please turn off the device immediately and contact our customer service.

11)The intensity level of the device will be initially set to Level 1. Adjust the energy level from low to high level. Please follow the recommended level protocol provided in the manual to choose the level that is suitable for your skin. After treatment, it is normal to have skin redness, suggest to use ice pad after treatment;

12)It is recommended to clean the treatment window after every session by using a clean cloth or cotton swab to gently wipe. This can avoid hair and sebum of carbon sticking glass and affecting energy of shooting light;

13)If necessary, only use a dry clean cloth to gently wipe the outer shell of the device. Never immerse the device or any of its parts in water and do not rinse it under the tap .

14)Do not use the device in wet surroundings (such as near bathtubs, showers, basins or other vessels containing water). If the device becomes wet, it should not be used;

15)Please store the device and its accessories in ventilated and dry environment and keep them dry and far away from pets, pests or children and lint, dust, light (including sunlight), etc.

16)Please do not put metal object into heat emission hole. Otherwise, there will be risk of electric shock;

17)In case the system is soaked please send it back to the factory for examination to determine whether it is repairable. In case the system is electrified when it is soaked, please power off before touching the system for your safety;

18)DO NOT disassemble, repair, or transform the device without authorization. If device cannot be operated normally, please send the system and accessories back manufacturer. Our qualified maintenance personnel will overhaul it.

19)This equipment is not needing special skills, training, and knowledge required of operator or responsible organization and environmental restrictions on locations of use.

20)People with reduced physical, sensory or mental capabilities or lack of experience and knowledge, may use the device if they have been given supervision or instruction concerning the safe use of the device and understand the hazards involved.

21)Regularly inspect the treatment window for damage: Do not use if the glass filter is cracked or missing. Do not use if the built-in reflector is missing or not in correct place.

22)Precautions to be taken in the event of changes in the performance of the device.

23)Precautions to be taken regarding the exposure of the device to reasonably foreseeable environmental conditions (e.g. to external electrical influences, electrostatic discharge, pressure or variations in pressure, acceleration, thermal ignition sources).

24) Instructions for installation, maintenance, pre-use checks and safe use, including clear warnings concerning precautions to avoid possible exposure to HAZARDOUS radiation or risk of fire.

25) When the device emits UV-A, UV-B or UV-C, care about the risks that PATIENTS will suffer from future potential carcinogenic occurrences.

26) To minimize visual disturbances to personnel, such as increasing the level of ambient light, applying spot lighting of the target area and use of protective eyewear.

27) Eye safety:

IPL Hair Removal Device emits flashes of Intense Pulsed Light. Through the Device has been tested and passed the international standards to ensure that its output is safe for your eyes, direct exposure is still potentially harmful to your eyes. The Device can only be flashed (light emitted) if both skin sensor (inside the treatment window) detect a valid skin reading. However:

--DO NOT look directly into the treatment window of the Device when the Device has electrical power and is switched on.

--DO NOT attempt to flash (emit light) the Device towards the eyes.

--DO NOT treat area around the eyes (eyebrows or eyelashes)

--ONLY attempt to flash the Device when the treatment window is in good contact with the area you wish to treat.

--WE RECOMMEND you look away from the Device when firing into your skin.

--WE RECOMMEND that you use protective glasses that comply with ISO 12609 when operating the device.

28) The consumption during normal use is electricity. For minimizing environmental impact during normal use, do not turn on when abnormal use, not to waste electricity.

29) The emissions during normal use are heat emission, the unused light energy emission and minor acoustic energy emission.

The fuse interrupting rating of switching power supply is 50A. The fuses of the switching power supply are used in the fault current is 50A, if the circuit interrupt current of the equipment is less than the fault current, the interrupt current is not used according to the fault current 50A.

1.3 Contraindications

1)DO NOT use if you are under the age of 18. The effects of using the system in this age group are unknown.

2)DO NOT use if the skin, in the area you want to treat, is tanned or burned. Use of the device on tanned or burned skin may result in temporary or even permanent whitening or darkening for the skin.

3)DO NOT use if you are pregnant, trying to become pregnant (conceiving) or breast-feeding. Hormonal changes occur during pregnancy or breast-feeding and could increase the risk of skin injury or skin sensitivity.

4)DO NOT use if you have epilepsy or light-induced seizures or migraines. The energy flashes might cause severe headaches and migraines.

5)DO NOT use on tattoos, dark brown or black spots (such as dark freckles, birthmarks, moles or warts). Tattoos or dark spots may absorb too much light energy which could cause skin swelling, temporary or permanent blisters or discoloration of the skin or whitening of the skin

6)DO NOT use on the head, ears, necks, nipples, genitalia, or around the anus. Use in these areas could cause skin swelling, blisters, or temporary skin discoloration or whitening of the skin.

7)DO NOT use if you have any active skin condition in the treatment area, including acne, sores, psoriasis, eczema, burns, herpes simplex, wounds or a rash. These conditions can make the skin more sensitive to the light treatment and could make the treatment uncomfortable or cause skin swelling, blisters or infection. If your skin condition has healed enough to use the device, then please consult with your doctor.

8)DO NOT use the device if you have a history of keloidal scar formation, a known sensitivity to light (photosensitivity) or are taking medication that makes the skin more sensitive to light, including non-steroidal anti-inflammatory agents, (e.g., aspirins, ibuprofens, acetaminophen), tetracyclines, phenothiazines, thiazide, diuretics, sulfonyleureas, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha- Hydroxi Acids(AHAS), Beta-Hydroxi Acids(BHAS),Retin-A, Accutane and/or topical retinoids.

9)DO NOT use if you have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc. These conditions make the skin more sensitive to the light treatment and could make the treatment uncomfortable or cause skin swelling, blisters or infection.

10)The IPL Hair Removal Device is MR Unsafe. Do not enter an MRI scanner room or an MR system. Doing so may result in severe injury or death.

1.4 EMC Statements

- IPL Hair Removal Device conforms to EMC requirements.
- User shall install and use in accordance with EMC information described in attached materials.
- During the use of the product, some of these RF emitters (e.g. RFID) in the intended environment of use might be concealed and give a potential harm to users. So it needs to pay attentions to the followings:
 - The portable RF communications equipment including antennas can effect medical electrical equipment, and it is used to no closer than 30cm (12 inches) to any part of the IPL Hair Removal Device;
 - Do not use the device in the present of or close to several common RF emitters (RF emitters that are known sources of electromagnetic disturbance such as diathermy, electrocautery, and RFIDs, walk-through metal detectors, and security systems (e.g., electromagnetic anti-theft systems, and metal detectors), these RF emitters have not been evaluated for electromagnetic compatibility with the subject device and it may be a source of electromagnetic disturbance and affect the safety of the device.
 - Do not apply the device in the presence of electronic monitoring equipment (e.g. cardiac monitor, ECG alarms), which may not operate properly when the electrical stimulation device in use.
 - The portable and mobile RF communication equipment may affect the performance of the device, so please avoid strong electromagnetic interference when using, such as do not use it near phone, microwave and so on.
 - It is prohibited to use this device in environments with specific medical EM radiation sources (such as electrocautery, MRI, electrosurgical units, and diathermy equipment). When using this electrical stimulation device, these electronic monitoring devices may not function properly, which could even lead to surgical delays or injuries. Therefore, it is prohibited to use the device in hospital operating rooms and clinical laboratories.
 - Do not use this device in environments with common EM radiation sources, including but not limited to radio frequency identification (RFID) readers, electronic security systems (metal detectors, electronic artificial surveillance (EAS)), wireless power transfer (WPT), and 5G cell signals. This may affect the performance of the device. Therefore, the device should be used in an enclosed room, away from the above common EM radiation sources. For example, do not use it near 5G transmitting and receiving devices (5G mobile phones, routers) or high-voltage transmission lines.
- Manual and manufacturer' s statement details are in the end of the article (Annex I).

The IPL Hair Removal Device is powered by Intense Pulsed Light (IPL) derived from professional IPL technology and it has the cooling function that effectively cools the skin. The IPL Hair Removal Device is an over-the-counter device indicated for the removal of unwanted hair. The device is also indicated for the permanent reduction in hair regrowth, defined as the long term, stable reduction in the amount of hair regrowing when measured at 6, 9 and 12 months after the completion of a treatment regime.

2.1 Packing List


When you open the IPL Hair Removal Device's gift box, you'll find it contains the following components:

IPL Main Device*1

AC Power Cord*1

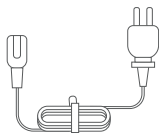
Storage Bag*1

User Manual*1

Goggles*1

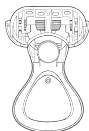
Razor*1


2.2 Product Components Diagram

EN

Here is the diagram of T32A:

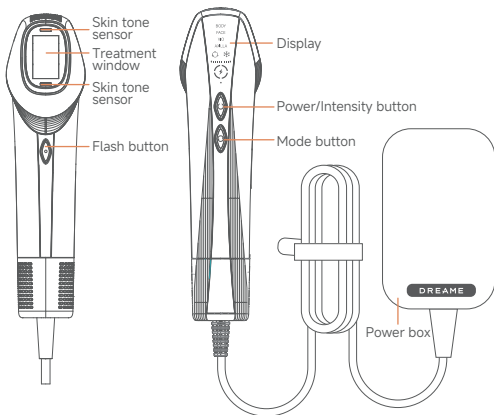- BODY — BODY AREA (One press 4 pulses)
- FACE — FACE AREA (One press 1 pulses)
- VIO — VIO AREA (One press 1 pulses)
- AXILLA — AXILLA AREA (One press 4 pulses)

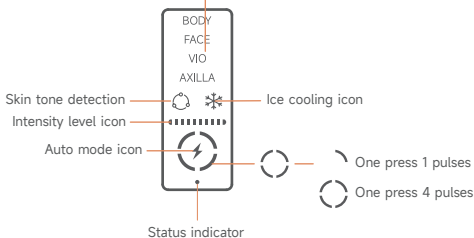IPL Hair Removal Device's Activation Buttons and Status Indicators:

1)Power/Intensity button:

· In either manual flash mode or automatic flash mode, short press the power button(for level adjustment) to switch from the auto skin tone detection level to level 1. A short press will cycle through the levels continuously in the following sequence: 1→2→3→4→5→6→7→8→9→10→Auto Skin Tone Detection level (Active).

· Pressing the button for over two seconds again, it will turn off the device.

2)Mode button:

· Short press the Mode button to switch to the treatment area (BODY/-FACE/VIO/AXILLA). The corresponding light emission count will be automatically matched after selecting a different area.

3)Flash Button:

· In manual flash mode, press this button to activate the flash lamp when the treatment window is in full contact with the skin. A short press triggers a single pulse or four pulses . Long press the flash button for 2 seconds to switch to Continuous Flash Mode. When entering the Continuous Flash Mode, flashes will be emitted at a certain interval, so please keep the sapphire treatment window in full contact with skin and keep moving the device. Under the Continuous Flash Mode, press the flash button for 2 seconds again to switch back to the manual flash mode.

· The area of FACE and VIO, a short press triggers a burst of 1 pulse.

· The area of BODY and AXILLA, a short press triggers a burst of 4 pulses.

Due to the high energy of body and axilla, When the device's skin tone detection function recognizes that it is being used on darker skin, the flash will not activate.

4)Status indicator light on the display:

| Status indicator | Mode | Meaning |
|------------------|--------------------|--|
| No Light | OFF | Device not working (Not connected to the AC Power cord). |
| Slow "flashing" | Activation/Standby | Device is in standby state (Device is activated after connected to the power supply). |
| "Steady" Light | Hair removal | Device can deliver a pulse (You can conduct "Hair Removal") due to only the full skin contact. |

3. How It Works

3.1 How Does IPL Work

IPL Hair Removal Device uses Intense Pulsed Light (IPL) technology to emit a specific wavelength of the light and delivers to the skin. The device is designed to help break the cycle of hair regrowth. Light energy is transferred through the skin's surface and is absorbed by melanin present in the hair shaft. The absorbed light energy is converted to heat energy (below the surface of the skin), which disables the hair follicle preventing further growth, so as to achieve effective hair removal.

3.2 What To Expect

Immediately Post-Treatment

You should not see any significant side effects from treatment (Refer to section 5.2 for more information on side effects).

In the first few weeks following initial treatments, you will still see some hair growing. These are likely to be hairs that were missed during treatment i.e. they were lying dormant and not in their growth phase (Anagen) during treatment, when IPL is most effective.

6 weeks into the 12 week plan

You should see a reduction in hair growth. However, many hairs may still not have been treated in their growth phase. It is important to continue the weekly treatments.

After the 12 week plan

Immediately after you have completed the 12 weeks plan you should see a significant reduction in hairs within the treated area. Any remaining hairs should be finer and lighter in colour. Continued monthly treatments, or treatments as required, should keep your skin hair-free.

4. How To Use

Preparation

IPL Hair Removal Device should ideally to be used on the small areas such as underarm and facial hair below the chin line, large area such as the legs.

4.1 Device Set-up

EN


1) Remove all components from the gift box;

2) Make sure the device is clean, dry and intact before beginning the treatment: You can clean the treatment window with clean cloth or cotton swab. This is done in order to remove hairs and dust from the treatment window;

3) Before each treatment, inspect the treatment window, AC Adapter cord for damage;

4) Plug AC adapter end into an electrical outlet and the cord into the device socket, the status indicator on the display will be blinking slowly;

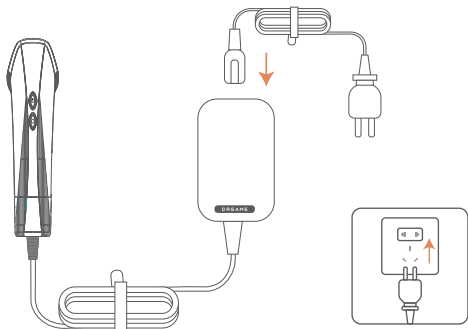IPL Hair Removal Device is now ready to be used.


NOTE:

- Place the IPL Hair Removal Device on a table or area near an electrical wall outlet. Move away any liquids from the area.
- In order to minimize the environmental impact during the expected service life, please throw the excess packaging materials into the dustbin after installing the equipment.

4.2 Preparing Yourself

1) On the basis of Fitzpatrick Photo typing Scale, identify your skin tone whether is suitable for; the details please refer to "4.4 Recommended Treatment Protocol" in this article.


2) Prior to use, make sure your skin is in good condition without any injury and wound; as well as your skin is shaved, clean and dry (free of any residue such as powders, creams, makeup, antiperspirants, or deodorants);

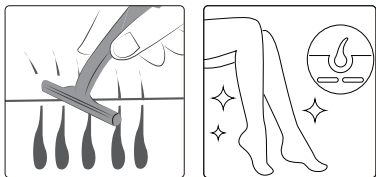4.3 Skin Sensitivity Test

If you are using the device for the first time or on a new body area, choose a small area close to the area you intend to treat, for the purpose of performing a skin sensitivity test.

- Use the LOW treatment level. Follow the treatment instructions below but make sure to emit only 2-3 pulses.
- After the skin test, wait 24 hours and check your skin for any reaction.
- If after 24 hours the skin appears normal and there is no redness or other skin reaction, then you are ready to perform a full treatment.

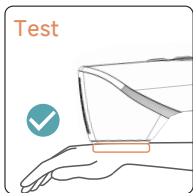4.4 Recommended Treatment Protocol

Applications:

Lip hair, armpit hair, body hair and legs, hair on places that affecting looks such as hair line on forehead and bikini area.

| Skin color Hair color | White | Ivory | Beige | Light Brown | Brown | Dark Brown |
|--|-------|-------|-------|-------------|-------|------------|
|  Black | ✓ | ✓ | ✓ | ✓ | ✓ | ✗ |
|  Dark brown | ✓ | ✓ | ✓ | ✓ | ✓ | ✗ |
|  Brown | ✓ | ✓ | ✓ | ✓ | ✓ | ✗ |
|  Golden | ✓ | ✓ | ✓ | ✓ | ✓ | ✗ |
|  Red | ✗ | ✗ | ✗ | ✗ | ✗ | ✗ |
|  White | ✗ | ✗ | ✗ | ✗ | ✗ | ✗ |
|  Gray | ✗ | ✗ | ✗ | ✗ | ✗ | ✗ |

✓ : Indicates applicable

✗ Indicates not applicable

NOTE:

- Choose the comfortable level according to your skin tolerance and increase level from one to ten levels gradually. If there isn't unsuitable or reaction, you could increase level gradually.

- The color of your skin may change, for example because your skin has been tanned: it may then be darker than the skin tone during the previous treatment. Therefore, it is important that you check the settings to be used before each treatment and set the settings according to different treatment areas. Some areas may be darker than others.

4.5 Usage Instructions


EN

1) Power On

Press and hold the Power/Intensity button for more than 2 seconds to turn on the device. After powering on, the device defaults to:

- BODY mode
- Automatic Skin Tone Recognition
- Ice Cooling function ON

The cooling function activates automatically and remains on throughout the treatment process.

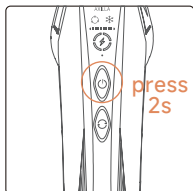
2) Select Mode

Short press the Mode button to switch between treatment modes:

- BODY: Suitable for large body areas such as arms and legs
- FACE: Suitable for facial areas, upper lip, and fingers
- VIO: Suitable for sensitive areas such as the bikini line
- AXILLA: Suitable for underarms and areas with coarser hair

The device will automatically match the recommended flash settings for each selected mode.

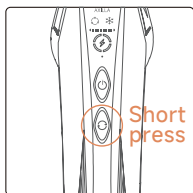3) Skin Cooling Preparation


After powering on, the Treatment window begins automatic cooling.

It is recommended to wait a few seconds until a noticeable cooling sensation is felt on the skin before starting treatment for enhanced comfort.

4) Eye Protection

Do not look directly at the flash light during use.
Always wear the goggles throughout the treatment.

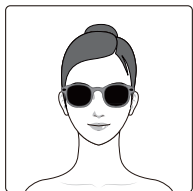5) Energy Level Control

Automatic Skin Tone Recognition

This mode is enabled by default after powering on

When the Treatment window is fully in contact with the skin, the device automatically detects skin tone and recommends an appropriate energy level (1–10)

Manual Level Adjustment

If manual adjustment is required:

In the current mode, briefly press the Power/Intensity button once to switch from Auto mode to Manual mode (Level 1–10 will be displayed)

Continue pressing the button to cycle through levels:

Auto → 1 → 2 → 3 → ... → 10 → Auto

Start from a lower level and gradually increase based on skin tolerance

To return to Automatic Skin Tone Recognition, continue cycling until Auto mode is selected

6) Flash Operation

EN



Manual Flash Mode

Keep the Treatment window fully in contact with the skin and press the Flash button once.

Continuous Flash Mode

For large areas, press and hold the Flash button for 2 seconds to activate continuous flashing, and glide slowly over the skin.

The ⚡ lights up, indicating continuous flash mode is active.

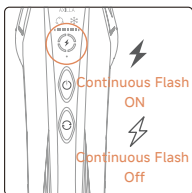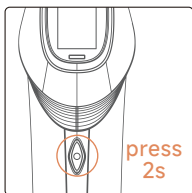7) During Treatment



Hold the treatment window at a 90 ° angle against the skin. When the window is fully in contact, the indicator light stays on, showing the skin is detected and the device is ready to emit energy.

After each emission, move the position of a light outlet in turn until the entire hair removal area is operated. Individuals can rest for a while according to different energy levels and the skin heat resistance. After the skin cools down, repeat the above steps for one or two times on the depilated area, but please do not repeatedly irradiate the same skin in a short period of time to avoid burning the skin.

Note: Duration of use: When using a single continuous treatment, avoid repeated use for more than three times on the same part.

Note: When the light outlet is moved away from the skin, no flash will be emitted.


Note: when in auto-flash mode, the user should move the device across the skin fast enough to avoid overlapping treatment areas. Not doing so could lead to overlapping treatment areas, and result in patient burns.

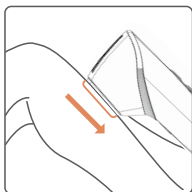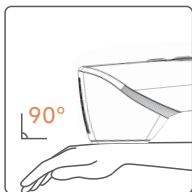8)After Treatment

Power Off


Press and hold the Power/Intensity button for more than 2 seconds to turn off the device. The fan will stop automatically after a short delay.

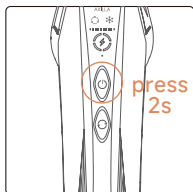Cleaning and Storage

After the device cools down:

- Clean the Treatment window according to the cleaning instructions
- Ensure no residual hair or dust remains
- Store the device and accessories in the packaging box
- Keep in a cool and dry place

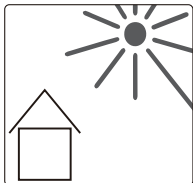4.6 Service life of the device

The valid maximum pulse counter of the device is up to 500,000 times. When the remaining pulse counter of the device shows '000000', it means that the device has reached the end of its service life. In this case, it needs to be replaced with a new one, and you can go to your local provider or contact the manufacturer directly.

4.7 ERROR reminder function

| Error Prompt Code | Defect Description |
|---------------------|-------------------------|
| ●○○○○○○○○● Flashing | Fan error |
| ●■■■■■■■■■ Flashing | low-voltage abnormality |
| ●■■■■■■■■■ Flashing | Temperature warning |

Solution: When the above error reminder appears, please turn off the device and disconnect the power, and then restart the device.

5.1 Possible Side Effects

When using IPL Hair Removal Device according to the instructions for use, side effects are rare. However, every aesthetic procedure including the usage of home use devices may cause some degree of side effects. If any of the following conditions (refer to section 5.2 Side Effect Table) appear and last more than 24 hours, you should consult your physician.

5.2 Side Effect Table

It is possible that the following side effects have emergence during use.

| Degree of Side Effects | Side Effects |
|----------------------------|---|
| Minor Side-Effects | Irritation / Redness / Sun Sensitivity |
| Rare Side-Effects | Burns / Scarring / Pigment Changes / Excessive Redness / Swelling |
| Extremely Rare Side-Effect | Infection / Bruising |

5.3 Cleaning, maintenance & storage

After each treatment, inspect the treatment window for damage. Do not touch or clean the treatment window immediately after use to allow it to cool down.

Cleaning

- Make sure to turn off and unplug main device before cleaning.
- Wipe the treatment window with a dry, lint-free cloth if any marks or black spots are visible.
- Never use acetone or other abrasive fluids, because they will damage your device.
- Keep IPL Hair Removal Device dry. Never submerge the device or any of its parts in water as if may lead to an electrocution risk.

Inspection

- After wiping, visual inspection the device. If the black spots remains, please wipe several times until no visual contamination.

Maintenance & Storage

- To avoid breakage or missing device and its parts, we recommend you store the device in its original box after each use.
- To avoid degradation of device performance, turn on the device at least once a year.
- Clean the device before you store it and store it in a safe and protective place.
- In order to minimize the environmental impact during the expected service life, please take good care of your equipment to avoid unnecessary damage. If it is damaged due to your own reasons, please place it in the e-waste recycling place reasonably, and do not discard it at will.

5.4 Trouble Shooting

EN

If you experience any problem with IPL Hair Removal Device, please refer to below information to try to resolve it:

| Problem | Solution |
|--|---|
| If IPL Hair Removal Device does not turn on: | Make sure that the power supply is plugged properly into the device and the electrical outlet. |
| If IPL Hair Removal Device does not trigger a pulse: | <ol style="list-style-type: none">1. Verify there is full contact with the skin. The flash button will activate a flash pulse only when the device is in full contact with the skin.2. In the case of a safety issue, the device will not trigger a pulse. |
| If you feel heat and hurt during the use of IPL Hair Removal Device. | <ol style="list-style-type: none">1. Shave and clean the hair in the expected treatment area before use.2. Choose one comfortable treatment level for yourself.3. DO NOT repeat treat the same area at the same time. |

If you still have problems with your IPL Hair Removal Device, please contact customer service.

6. Information

6.1 Technical Specifications

| Basic Specifications | |
|--|---------------------------------------|
| Product Name | IPL Hair Removal Device |
| Model | T32A |
| Trade Mark | DREAME |
| Level of Protection | BF Type, Class II |
| Service life of the device | 500,000 times (the maximum shot time) |
| Shelf life of the device | 5 years |
| Method(s) of sterilization | N/A |
| Suitability for use an OXYGEN RICH ENVIRONMENT | N/A |
| Mode of operation | Continuous operation |

| Output Specifications | |
|--|--------------------------------|
| Type of Light | Intense Pulsed Light (IPL) |
| Energy Medium | Xenon Arc flashlamp |
| Light Wavelength Range | 510~1200nm |
| Pulse Duration | 0.4ms~0.92ms Multipulse |
| Spot Size (Window size of IPL) | 5.4cm ² |
| Treatment Level | Level 1 to 10 |
| Energy output(J)±20% | 10.0~35.0J |
| Energy Density(J/cm ²)±20% | 1.85~6.48J/cm ² |
| Risk Group | Exempt |
| <p>Caution: Use of controls or adjustments or performance of procedures other than those specified herein may result in HAZARDOUS radiation exposure.</p> <p>The product fulfills the requirements of IEC60601-2-57:2011and IEC 60601-2-83:2019/AMD1:2022.</p> | |
| Power Specifications | |
| Input | AC 100~240V, 50/60Hz, 3.5A MAX |
| Usage Specifications | |
| Operation Temperature | 0°C~30°C (32°F~86°F) |
| Operation Humidity | 25%~75% |
| Atmospheric Pressure | 80kPa~106kPa |
| Storage Temperature | 0°C~50°C (32°F~122°F) |
| Storage Humidity | 25%~75% |
| Atmospheric Pressure | 80kPa~106kPa |

6.2 Disposal

DO NOT throw away the device with normal household waste at the end of its life, but hand it in to an official collection point for recycling (Contact your local town or city officials for recycling information). By doing this, you help to preserve the environment.

Note:

- Consumption during normal use just electric energy;
- Emissions during normal use are acoustic energy and heat;
- The location of components and parts within the device that contain stored energy near main board; To disassembler IPL Main Device, you need to place it more than 5 minutes to discharge the capacitor.

6.3 Disclaimer

Dreame Trading (Tianjin) Co., Ltd. reserves that the final explanation of IPL Hair Removal Device, any third party applies device related information without Dreame's authorization is not allowed, we will take legal action. In any case, Dreame Trading (Tianjin) Co., Ltd. is not responsible for any consequences caused by improper use of consumers.

6.4 Support

Dear Customers:

In order to provide you with more convenient services, when you buy IPL Hair Removal Device, please consult according to below contact information if you have problems.

Manufacturer:

Dreame Trading (Tianjin) Co., Ltd.

Address: Room 2112-1-1, South District, Finance and Trade Center, No.6975 Yazhou Road, Dongjiang Bonded Port Area, Tianjin Pilot Free Trade Zone, 300463 Tianjin, P.R. China

Customer Service:

For any questions or issues, please contact us at:

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Manufacturer: Dreame Trading (Tianjin) Co., Ltd.