Dreame Chrona Photon Anti-Aging Skincare Mask

User Manual

DREAME

EN	User Manual.....	001
RU	Руководство пользователя	015
DA	Brugsanvisning	029
PL	Podręcznik Użytkownika	043
FI	Käyttöohje.....	057
NO	Bruksanvisning	071
FR	Manuel d'utilisateur	085
GR	Εγχειρίδιο χρήστη	099
NL	Gebruikershandleiding	113
DE	Benutzerhandbuch	127
ES	Manual de usuario	141
IT	Manuale dell'utente	155
PT	Manual do utilizador	169
SV	Bruksanvisning	183
TR	Kullanım Kılavuzu	197
HE	שמתשמל גירדמ	211
AR	مدخست سمدل ليلد	225
ESLA	Manual do utilizador	239
TH	คู่มือผู้ใช้.....	253
VN	Hướng dẫn sử dụng	267
ID	Panduan Pengguna	281
BM	Panduan Pengguna	295
TW	使用者手冊	309

Warning

To prevent electric shock, fire, or other accidents, carefully read this user manual before use and retain it for future reference.

- Before using the Dreame Chrona Photon Anti-Aging skincare Mask you must carry out a sensitization test (see section light sensitivity test). DO NOT use the Dreame Chrona Photon Anti-Aging skincare Mask if you are sensitive to Dreame Chrona Photon Anti-Aging skincare Mask. This unit is intended to be used in the home healthcare environment, the patient can be an intended operator.
- DO NOT use the Dreame Chrona Photon Anti-Aging skincare Mask in any way that is not exactly as specified in these instructions for use. Use of the Dreame Chrona Photon Anti-Aging skincare Mask other than exactly as specified in these instructions for use has not been tested and the risk is unknown. This equipment is not intended for use by persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the equipment to avoid the risk of fire or burns.
- Turn OFF and remove the Dreame Chrona Photon Anti-Aging skincare Mask immediately if the device feels hot or if your skin becomes uncomfortably warm.
- This product is wireless during use. A charging cable is only connected when charging.
- During charging, the cable may present a strangulation hazard. Keep the device and charging cable out of reach of children at all times.
- DO NOT use the Dreame Chrona Photon Anti-Aging skincare Mask without first removing makeup BEFORE treatment. Some ingredients in these products may cause light sensitization or may reduce the effect of the treatment. If you use the Dreame Chrona Photon Anti-Aging skincare Mask without removing these products, you may cause a skin reaction, or you may reduce the effect of the treatment. This equipment is not intended for use by children. Children should be supervised to ensure that they do not play with the equipment to avoid the risk of fire and burns.
- NEVER use Dreame Chrona Photon Anti-Aging skincare Mask near water or where water can be spilled onto the equipment. DO NOT use the Dreame Chrona Photon Anti-Aging skincare Mask if you think it has come into contact with water. Do not subject the equipment to heavy shocks to avoid risk of damage to the lamp.
- Keep the Dreame Chrona Photon Anti-Aging skincare Mask away from HOT surfaces. If needed, contact the manufacturer for assistance in setting up, using, or maintaining the ME EQUIPMENT or reporting unexpected operation or events.

- EN**
- DO NOT use the Dreame Chrona Photon Anti-Aging skincare Mask if the mask, remote control, or the USB Type-C charging cable is damaged in any way. Do not leave the device unattended when it is powered on to avoid the risk of fire, burns, or other hazards.
 - When removing from minimum / maximum storage temperature, it is recommended to recover for 2 hours before use.
 - Do not look directly or shine light in the eyes while in operation.
 - DO NOT use the Dreame Chrona Photon Anti-Aging skincare Mask nearby the device with wireless technology.

Skincare Restrictions (Before Use)

To ensure safety & phototherapy efficacy, avoid facial products with these ingredients:

- Vitamin A & derivatives (retinol, retinaldehyde, retinyl esters, etc.)
- Acids (fruit acids, high-concentration salicylic acid, etc.)
- Irritants (high-concentration alcohol, complex fragrances, artificial colorants)

Recommended Steps

1. Cleanse face thoroughly; keep skin dry and residue-free before use
2. Post-treatment: Apply mild moisturizing/repair products (skin absorption is enhanced)

Age Eligibility

For adults 18 years old and above only

Contraindicated Groups

Do not use, or consult a doctor first, if you:

- Have open facial wounds, eczema, dermatitis, acute rosacea, severe acne or active skin conditions
- Have photosensitive diseases (e.g., lupus erythematosus, solar dermatitis) or take photosensitive medications (e.g., certain antibiotics, antihypertensives)
- Are pregnant, breastfeeding or under 18
- Have facial/head implants (pacemakers, stents, defibrillators, etc.)
- Underwent laser, chemical peel or microneedling within 1 month (wait for full skin recovery)
- Are unsure: Consult a doctor/dermatologist first

Usage Tips

- Ensure mask fits snugly to skin; avoid contact lenses
- Strengthen daily sun protection during use to enhance results and prevent pigmentation
- Stop use immediately if discomfort occurs; consult a dermatologist if symptoms persist

Disclaimer: This is a home-use beauty device. Not a substitute for professional medical diagnosis or treatment.


Precautions

- This device contains silicone, plastic components, and hook-and-loop fastener materials. Individuals with known allergies to any of these materials should discontinue use and consult a healthcare professional.
- This device is equipped with a memory function. It automatically recalls the last-used treatment mode, intensity level, and treatment time each time the device is powered on.
- The Dreame Chrona Photon Anti-Aging skincare Mask cannot be used in charging mode.
- Before each use, with the device switched off, clean the face mask.
- DO NOT use any other cleansing products other than those stated in the cleaning instructions. The use of other cleansing products may damage the Dreame Chrona Photon Anti-Aging skincare Mask.
- Store the Dreame Chrona Photon Anti-Aging skincare Mask at room temperature in dry conditions, away from direct sunlight.
- When removing from minimum / maximum storage temperature, it is recommended to recover for 2 hours before use.
- Do not look directly or shine light in the eyes while in operation.

Electrical safety

- Do NOT damage, modify, bend, twist, or heat the USB Type-C charging cable. A damaged charging cable may lead to fire, electric shock, or malfunction.
- DO NOT place or operate the device near water or liquids. It may cause fire or electric shock.
- This product uses a magnetic, detachable battery pack. Charge the battery pack fully before using the device for the first time.
- Only use the designated battery pack provided with the device. Do NOT attempt to open, modify, replace, or use any non-approved battery packs, as improper handling may lead to overheating, leakage, or explosion.
- DO NOT dispose of the battery randomly, follow your local disposal procedure when disposing of it.
- Use only a DC 5V = 1A charging adapter to charge this device (maximum output power 3.5W). Do NOT use chargers of unknown origin, poor quality, or with an output power exceeding 30W, as this may cause overheating or damage to the product. For maximum safety, it is recommended to use a charging adapter that complies with IEC 60601-1:2005 + AMD2:2020, with the following rating: Input: 100-240VAC Output: DC 5V = 1.0A.

EN Product Overview

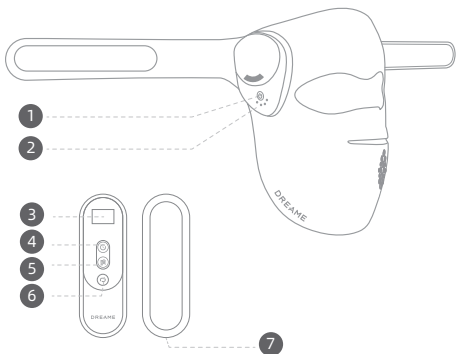- 1 Turning on/off
- 2 Battery Indicator:
 - 4 Green Lights: 100% Charge (Fully Charged)
 - 3 Green Lights: 75% Remaining
 - 2 Green Lights: 50% Remaining
 - 1 Green Light: 25% Remaining
 - Flashing Red Light: <10% Low Battery
- 3 LCD for checking the selection mode/Intensity /Timer
- 4 Timer setting button
5 mins/10 mins/ 15mins
- 5 Intensity Levels button:
 - Level I (Low)
 - Level II (Medium)
 - Level III (High)
- 6 On/Off & Treatment modes selection button:
 - Mode 1: Rejuvenate
 - Mode 2: Brighten
 - Mode 3: Balance
 - Mode 4: Calm
 - Mode 5: Restore
- 7 Type-C port

The product body is marked with "Do not use while bathing or showering".

All product images, accessories, and user interface illustrations in this manual are for reference only. Due to product updates and upgrades, the actual product may differ slightly from the illustrations. Please refer to the actual product.

Getting started

Device Set-up

Getting started with your Dreame Chrona Photon Anti-Aging skincare Mask is simple.

1. Unpack the device and all included accessories.
2. This product uses two magnetic, detachable battery packs (each 1000mAh, total 2000mAh). Before first use, fully charge both battery packs and the remote control. Refer to Charging Instructions for details.

Light sensitivity test


Follow these instructions to see if you are sensitive to light:

1. Press the Power button on the LED mask to start the device. Use the remote control to switch the device On/Off and to choose between the available treatment mode(color).
2. Place the device on a section of the forearm, maintaining contact for the entire three-minute treatment cycle.
3. Press the "Treatment modes selection button" again to choose a different treatment mode(color) for another three-minute treatment on the same forearm area.
4. Wait 24 hours. You should NOT experience any adverse reaction during the next 24 hours.
5. An adverse reaction is a pain or discomfort in the treatment areas or any pro- longed skin redness.

Charging Instructions (Battery Packs & Remote Control)

This device is equipped with two detachable 1000mAh battery packs and a wireless remote control (280mAh).

Please follow the instructions below before first use:

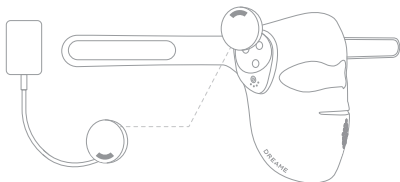Charging the Battery Packs

- Remove the battery pack from the mask before charging.
- Each battery pack requires up to 2 hours to reach a full charge.
- Charging time may vary depending on the remaining battery level.

Charging the Remote Control

- The wireless remote control includes a 280mAh rechargeable battery.
- The maximum charging time is up to 2 hours.
- Charging time may vary depending on the remaining battery level.

Lithium Battery Maintenance Tips:


- Avoid fully draining the battery before recharging.
- Avoid overcharging the battery and try to limit the charging time to a maximum of 8 hours.
- If the device is not in use for an extended period, it is recommended to charge the battery at least once every three months.

How to use

Preparation

1. Wash your face clean. (If you make up, please cleanse by cleansing foam.)
2. It can be used day or night. You can fit it into your current skincare routine after cleaning and/or exfoliating your skin.
3. Be sure to fully charge the battery before use.

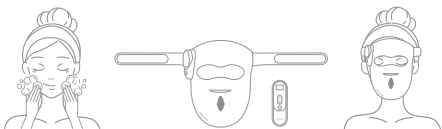1. Shortly press the button of the LED mask to turn on the power. (Available when blinking green)
2. After adjusting the hair band moderately, wear it so that the inner side of the LED mask is on the face.

Tips for adjusting the headband:

1. When sitting or standing, the correct state should be minimal pressure on the eye socket, allowing the silicone eye pad to fit the eye socket, with the headband fixed the product.
2. If you feel the silicone eye pad is exerting too much pressure on the eye socket, loosen the headband on both sides.
3. If the silicone eye pad does not fit the eye area and light can enter through the gap between the eye pad and eye socket, tighten the headband on both sides.
4. If you feel the mask is sliding down and cannot maintain the correct position on your face, tighten the headband.
5. Select the desired mode using the buttons on the remote control.
6. If you feel uncomfortable while wearing the device, you can try to adjust the strap or pick it off until you feel better.
7. The product will be turned off automatically when the set mode and the corresponding time finished.

Cautions for Use

- Move the remote control within 2 m of the LED mask.
- Please do not act violently when wearing the product.
- When wearing the product, relax in a comfortable position on the couch or bed.

Operation instruction

Remote control

1. Connect the Type-C port of the remote control using the supplied charging cable.
2. Plug the USB end of the cable into a 5V 1A power adapter to charge.
3. Flashing "88" Display: Charging in progress
Solid "88" Display: Fully charged
Unplug the charging cable before use.
4. When the remote control battery is low, the "88" display will flash until the screen turns off.

Since the LED mask uses a radio frequency remote control (the remote emits radio frequency signals, which the main unit receives), users need to pair the remote control with the LED mask when using it for the first time. The pairing steps are as follows:

1. Press and hold the power button on the LED mask for 5 seconds to enter pairing mode. The mask's indicator light will switch from green to red and then flash rapidly between green and red.
2. Press and hold the remote's power button for 3 seconds to turn on. A continuous beep and steady green indicator light confirm successful pairing.

Operating Instructions

Treatment Mode Selection

Press the ON/OFF & Treatment Modes Selection Button (which includes both the power button and treatment mode selection function) to choose from the following treatment modes:

Mode 1: Rejuvenate Red Light 630nm + Near Infrared Light 850nm

Benefits:

Targets signs of skin aging by stimulating collagen production, enhancing skin elasticity, and promoting cellular repair. This mode helps restore a smoother, firmer, and more radiant complexion, leaving the skin looking rejuvenated and youthful.

Mode 2: Brighten Yellow Light 590nm

Benefits:

Soothes and revitalizes the skin, providing exceptional recovery benefits after sun exposure or fatigue. It helps brighten and restore a vibrant, even-toned complexion, promoting a refreshed and youthful glow.

Mode 3: Balance Green Light 530nm

Benefits:

Calms and balances sensitive or stressed skin, reducing redness and promoting an even skin tone. This mode provides a naturally balanced and refreshed appearance, ideal for those with reactive or sensitive skin.

Mode 4: Calm Blue Light 465nm

Benefits:

Often used in beauty and skincare routines, this mode helps soothe the appearance of visible redness, reduce the look of oiliness, and support a clearer-looking complexion. It helps keep the skin feeling fresh and balanced, contributing to healthier and more refined-looking skin.

Mode 5: Restore Full-Spectrum Mode: Yellow 590nm + Green 530nm + Red 630nm + NIR 850nm

Benefits:

This multi-light therapy mode offers a comprehensive approach to a variety of skincare concerns, such as anti-aging, skin brightening, pigmentation reduction, and overall skin clarification. The combination of different wavelengths provides a holistic treatment to rejuvenate the skin, even out skin tone, fade dark spots, and enhance skin texture. Ideal for daily, all-in-one facial care, it delivers both restorative and preventive benefits for a youthful, glowing complexion.

Intensity Levels

Press the Intensity Levels Button to select the intensity:

Level I (Low)

Level II (Medium)

Level III (High)

Timer Settings

Press the Timer Setting Button to select the timer:

· 5 mins

· 10 mins

· 15 mins

Memory Function

Upon first use, the mask will default to **Mode 1: Rejuvenate, Intensity Level III(High), and Timer 10 mins.**

After each use, the mask will remember the last used treatment mode, intensity, and timer settings.

Notes

Feature

Description

Power On/Off

Press and hold the ON/OFF button for 5 seconds to enter standby mode (indicator light turns green).

Treatment Mode Switching

In standby mode, short-press the remote control button to switch between treatment modes: Rejuvenate → Brighten → Balance → Calm → Restore.

Intensity Settings

In standby mode, short-press the remote control button to switch intensity: Low → Medium → High.

Timer Settings

In standby mode, short-press the remote control button to select timer: 5mins → 10mins → 15 mins.

Auto Shutdown

1. In standby mode, the mask will automatically turn off after 10 minutes of inactivity.
2. The mask can be turned off by pressing the button on the mask.
3. The mask will automatically turn off once the treatment duration ends.

Memory Function

The mask will remember the last used treatment mode, intensity, and timer settings, and will restore these settings when powered on again

Additional Instructions

1. When the set time ends, the mask will automatically turn off.
2. When the remote control is turned off, the mask will also shut down.
3. After selecting a mode with the remote, do not use the mask buttons, as this resets the timer to 10 minutes and causes a display mismatch. Start treatment directly after selecting the mode with the remote control.

Care and maintenance

1. Store the Dreame Chrona Photon Anti-Aging skincare Mask at room temperature in dry conditions. Use the carry bag provided.
2. Avoid storing the Dreame Chrona Photon Anti-Aging skincare Mask in dusty conditions, in direct sunlight or in damp and humid conditions.
3. Dreame Chrona Photon Anti-Aging skincare Mask contains no serviceable parts.
4. Before each use, with the device switched off, clean the face mask following the cleaning instructions below.
5. Avoid bending or twisting the silicone mask excessively.
6. Do not apply cosmetics directly onto the mask surface.
7. Keep the product away from oils, lotions, and perfumes.

Cleaning the Dreame Chrona Photon Anti-Aging skincare Mask

1. The mask body is made of silicone. After each use, gently wipe the surface with a soft, lint-free cloth dampened with clean water.
2. If needed, clean with mild, alcohol-free, fragrance-free neutral soap diluted with water.

Disposal

DO NOT dispose any part of the Dreame Chrona Photon Anti-Aging skincare Mask at a domestic waste facility. Dreame Chrona Photon Anti-Aging skincare Mask and remote control contains a Lithium-ion polymer battery. For the protection of the environment, please package and return this unit to the manufacturer or their representative at the end of its working life or safely dispose of the device following your local guidelines.

Common malfunctions and solutions

Issue	Solution
<p>The screen on the remote does not light up, and the mask cannot be controlled.</p>	<p>Check if the remote control has power. If the battery is dead, please charge the remote using the included charging cable and a 5V 1A adapter. If the remote still does not work after charging, please contact us to arrange a replacement.</p>
<p>After turning on the mask, the remote screen lights up, but the mask cannot be controlled.</p>	<p>Please reset both the remote control and the mask. Pairing Instructions: a. Press and hold the mask's power button for 5 seconds. The indicator will first show a green light and then a flashing red light, signaling that the device has entered pairing mode. b. Simply press the ON/OFF button on the remote control to complete the pairing process.</p>
<p>The mask cannot be charged.</p>	<p>a. Replace the charging cable. If it is indeed a charging cable issue, you can contact us for a replacement of the charging cable. If the problem still cannot be solved, please contact us to replace the device. b. Check whether the charging adapter to which you connect the unit is live or whether the voltage indicated on the unit corresponds to the local mains voltage.</p>

EN Basic Parameters

Product Name	Dreame Chrona Photon Anti-Aging Skincare Mask
Model	AB01A
Battery specifications	Mask Battery Packs • Type: 3.7V / 1000mAh Lithium Polymer Battery (per battery pack) • Quantity: 2 detachable battery packs (total capacity: 2000mAh)
Product Dimensions	Mask : 634*233.8*5.5mm Remote Control: 133.6*42.3*19mm
Remote Control	Type: 3.7V / 280mAh Lithium Polymer Battery
Net Weight	Mask: About 338g Remote Control: About 72g
Adapter Input	5V $\overline{=}$ 1A
Operating Time	• One fully charged battery pack: Up to 60 minutes • Two battery packs (combined usage): Up to 120 minutes total
Charging time	• One battery pack: Up to 2 hours • Two battery packs (charged separately): Up to 4 hours total • Remote control: Up to 2 hours

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